


# Zirkonzahn®

*Human Zirconium Technology*

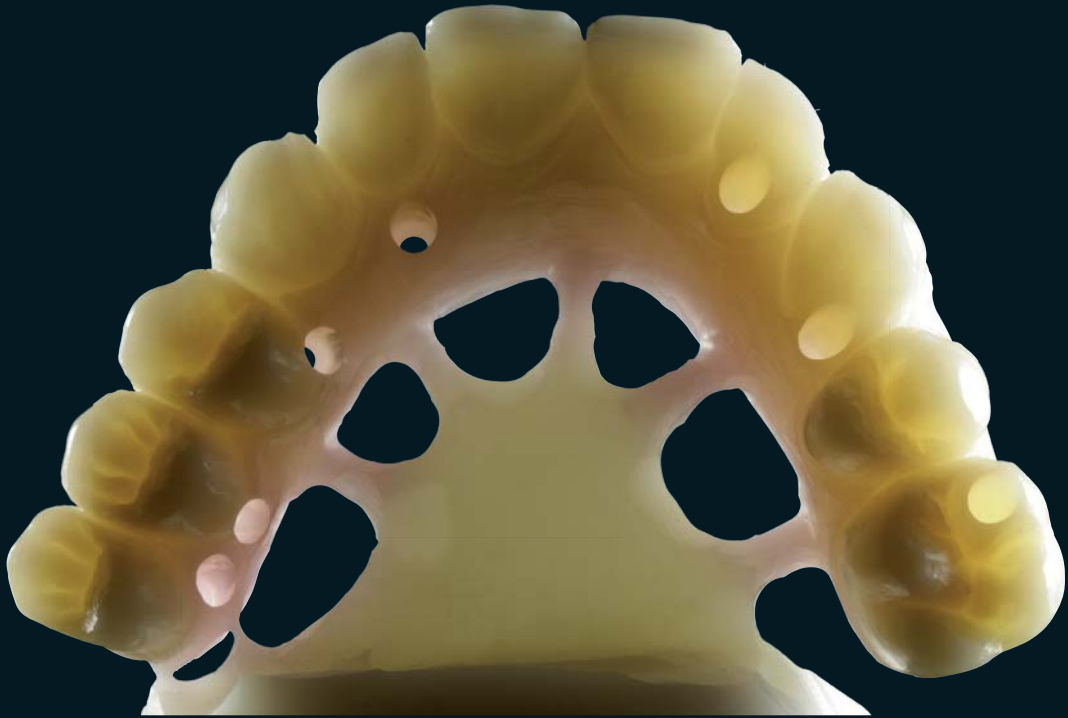

**100% PRETTAU ZIRCONIA**

## *Prettau Zirconia - 2 full zirconia bridges*

The colouring technique is a key factor for a three-dimensional colour effect and a natural result. By means of the Prettau colouring technique, at first the tooth area and then the gums are coloured. The Prettau intensive colours Incisal grey, Incisal violet and the Colour Liquid fluorescence were used in the anterior tooth region. A thin dentine layer was applied at first on the framework, before the ceramic layer. For a natural appearance of the gums, a wide as possible, varied colour spectrum should be used for the gingival layer. The anterior teeth were then simply veneered with a thin enamel layer. The gingival area should be layered at the same time as the labial area. Finally, a manual polish provides the finished work with a last natural touch.



## *The Prettau Bridge*

Jürgen Auffarth, Austria

In cooperation with ZTM Christian Hannker, Prosthetic dentistry laboratory Bellmann & Hannker (Germany)




# Zirkonzahn®



## HUMAN ZIRCONIUM TECHNOLOGY

ZIRKONZAHN World Wide - Tel +39 0474 066 680 - [info@zirkonzahn.com](mailto:info@zirkonzahn.com) - [www.zirkonzahn.com](http://www.zirkonzahn.com)

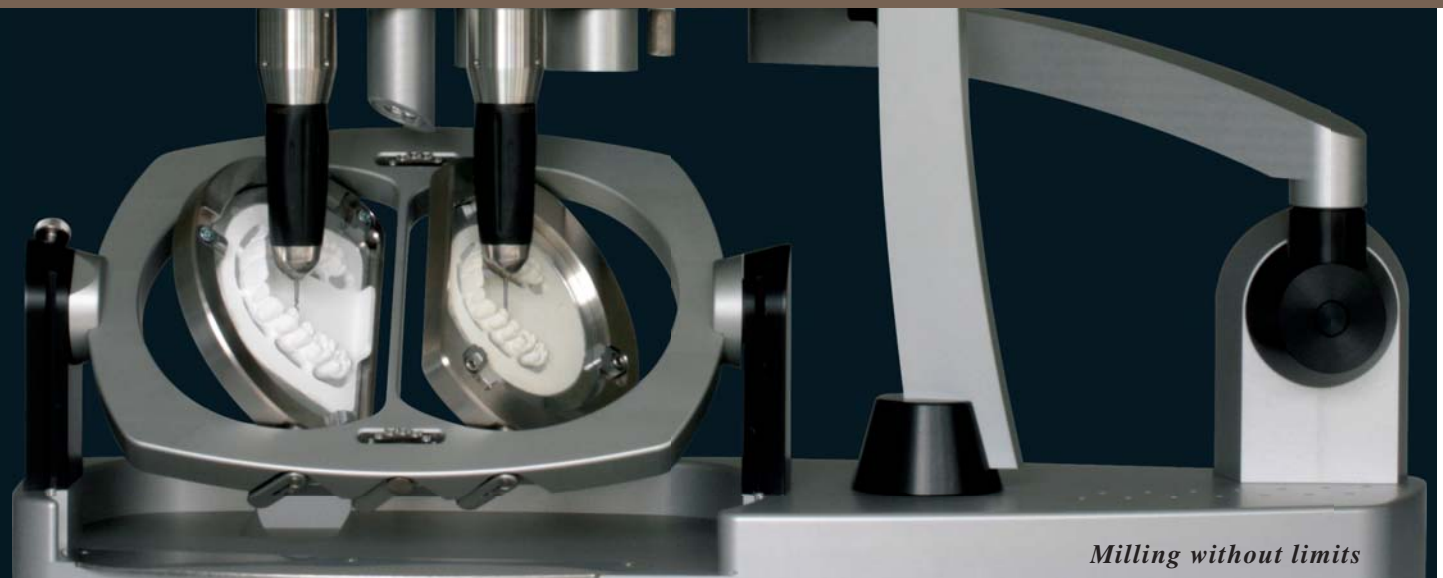

*Milling without limits*